

KANMED[®] BABYBED

KANMED BABY WARMING SYSTEM

Kanmed BabyBed offers you many choices. See how you can configure the bed to meet your requirements.



Water at the right temperature provides optimal warming

37°C



Configure your Kanmed° BabyBed



Water at the right temperature
provides optimal warming

The Kanmed° BabyBed is modern, practical and made for the care of newborn and premature babies in medium and low intensive care situations. The BabyBed was developed together with healthcare professionals to ensure a high level of ergonomics. Kanmed° BabyBed gives you unique push down side walls, ease of use and is very easy to clean. The Kanmed° BabyBed is manufactured in Sweden and comes with 5 years limited warranty.

The Kanmed° bed is often used with the Kanmed° BabyWarmer creating the Kanmed BabyWarming System. With the Kanmed° BabyWarming System you will need fewer incubators and reduce the overall costs. The Kanmed° BabyBed has a unique selection of options and accessories. By selecting models and accessories you will get the bed that fits your needs.

SELECT BED MODEL WITH OR WITHOUT RAIL

With surrounding plastic rail

- ☐ BB100 Electrical height adjustment
- ☐ BB300 Fixed height
- ☐ BB400 Electrical height adjustment
- ☐ BB400F Fixed height

Without surrounding rail

- ☐ BB101 Electrical height adjustment
- ☐ BB301 Fixed height
- ☐ BB401 Electrical height adjustment
- ☐ BB401F Fixed height

With surrounding metal rail

- ☐ BB102 Electrical height adjustment
- ☐ BB302 Fixed height
- ☐ BB402 Electrical height adjustment

TENT OPTIONS

Various outside and inside colours are available, see www.kanmed.se for the latest designs.

- ☐ BB160 (standard) (A)
- ☐ BB162 (B)
- ☐ BB161 (C)
- ☐ BB460 (twin standard) (D)



A



B




C



D

ACCESSORIES

- ☐ 800-0518 Drip Pole ø 25 mm, L 66 cm (A)

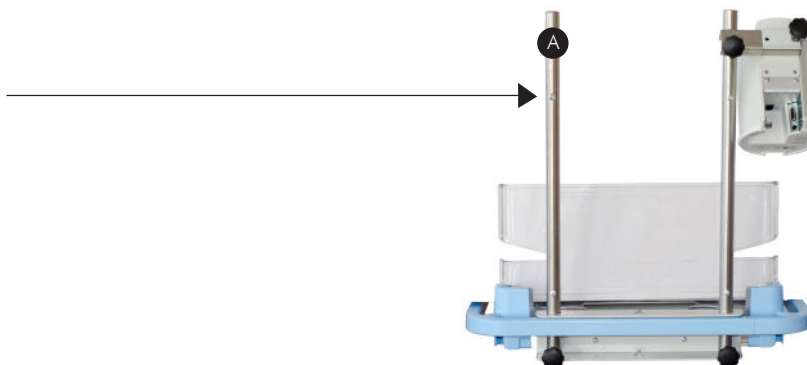

STANDARD BED DELIVERY

All Beds include: Tent BB160, (BB460 for Twin Bed), Tent Pole, Mattress (sealed), Side walls, one Storage Box BB170, Power cable and user manual.



On all Kanmed® beds

Twin version with
2 individually
adjusted Kanmed®
Baby Warmers
and foot pedal front



With rail



Without rail



ACCESSORIES


- BB171 Bottle Holder
- BB173 Trendeleburg kit for +- 17 degree tent
- 800-0502 Cable for BabyWarmer control unit (outlet on Pillar, one cable less on the floor)
- BB170 Storage box in front (extra)
- BB180 Shelf 33x33 cm edges 3cm (plastic coated steel)
- BB190 Shelf 33x33 cm edges 5cm (plastic coated steel)

ELECTRICAL HEIGHT ADJUSTMENT


- BB175 Foot pedal left (standard)
- BB175 Foot pedal right (note on your order)
- BB176 Foot pedal front (note on your order)

WHEELS

- 700-0874 Standard 125 mm (A)
- 700-0884 Double wheel 95 mm (B) (Lowers bed about 4 cm)

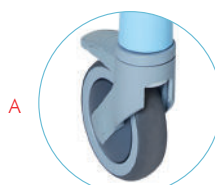

A

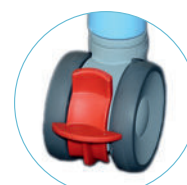

B

KANMED[®] BABYBED

	BB100 / 101 / 102 Standard BB300 / 301 / 302 Fixed height		BB400 / 401 Twin / 402 BB400F / 401F Fixed height Twin	
Outer dimensions	L 85 cm W 63 cm With rails	L 83 cm W 52 cm Without rails	L 95 cm W 80 cm With rails	L 93 cm W 70 cm Without rails
Wheel base / Wheels	87 cm x 64 cm / Diameter 12,5 cm alt. 9,5 cm, with brakes			
Bed Tray dimensions	L 69 cm W 46,5 cm H 9 cm		L 76,5 cm W 64 cm H 9 cm	
Height from floor to Bed Tray bottom	Min 78 cm, Max 108 cm. Fixed height is 80 cm (9,5 cm wheels lower all beds about 4 cm)			
Weight	About 46 kg Fixed height 36 kg		About 52 kg Fixed height 42 kg	
Tilt	About 17 degree (option ± 17 degree)			
Tent / Mattress	Washable 60°C (90°C) / Soft foam, sealed watertight cover			
Storage Box / Shelf	L 32 cm W 33 cm H 12 cm, Lid 32 x 33 cm / 33 x 33 cm			
Voltage BB100 / 400	230V 50Hz / 60Hz or 115V 60Hz, 150VA. Protection BB100 / 400: Class I			

► Note that the Baby Warmer may be used separately in any suitable bed.


KANMED[®] BABYWARMER

It is the soft warm
Water Mattress that makes
the Kanmed Baby Warming
System unique.



The Kanmed Baby Warming system is the most baby friendly way to care for your babies. Only a soft warm water mattress can substitute the warming efficiency of a modern incubator. When the baby is stable and the skin has matured it can be moved to the open Kanmed Baby Warming System. Parents can now easily get involved and practice the care of their new-born or premature baby.

The important parent baby bonding and breast feeding is now much simpler. The Kanmed Baby Warming System is used in hospitals worldwide as a safe, cost effective and baby friendly alternative to continued treatment in the incubator. Thereby the amount of incubators required can be reduced saving money for other investments in the Neonatal Unit. The system is well documented in published medical studies.

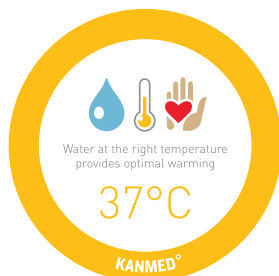

More than 10 000 units in use worldwide
that have warmed more than 2 million babies